

ACCEPTED MANUSCRIPT

## Chest electrical impedance tomography measures in neonatology and paediatrics – a survey on clinical usefulness

To cite this article before publication: Inez Frerichs *et al* 2019 *Physiol. Meas.* in press <https://doi.org/10.1088/1361-6579/ab1946>

### Manuscript version: Accepted Manuscript

Accepted Manuscript is “the version of the article accepted for publication including all changes made as a result of the peer review process, and which may also include the addition to the article by IOP Publishing of a header, an article ID, a cover sheet and/or an ‘Accepted Manuscript’ watermark, but excluding any other editing, typesetting or other changes made by IOP Publishing and/or its licensors”

This Accepted Manuscript is © 2018 Institute of Physics and Engineering in Medicine.

During the embargo period (the 12 month period from the publication of the Version of Record of this article), the Accepted Manuscript is fully protected by copyright and cannot be reused or reposted elsewhere.

As the Version of Record of this article is going to be / has been published on a subscription basis, this Accepted Manuscript is available for reuse under a CC BY-NC-ND 3.0 licence after the 12 month embargo period.

After the embargo period, everyone is permitted to use copy and redistribute this article for non-commercial purposes only, provided that they adhere to all the terms of the licence <https://creativecommons.org/licenses/by-nc-nd/3.0>

Although reasonable endeavours have been taken to obtain all necessary permissions from third parties to include their copyrighted content within this article, their full citation and copyright line may not be present in this Accepted Manuscript version. Before using any content from this article, please refer to the Version of Record on IOPscience once published for full citation and copyright details, as permissions will likely be required. All third party content is fully copyright protected, unless specifically stated otherwise in the figure caption in the Version of Record.

View the [article online](#) for updates and enhancements.

# Chest electrical impedance tomography measures in neonatology and paediatrics – a survey on clinical usefulness

Inéz Frerichs<sup>1</sup> and Tobias Becher<sup>1</sup>

<sup>1</sup> Department of Anaesthesiology and Intensive Care Medicine, University Medical Centre Schleswig-Holstein, Campus Kiel, Kiel, Germany

E-mail: frerichs@anaesthesie.uni-kiel.de

## Abstract.

*Objective:* Chest electrical impedance tomography (EIT) is currently applied in neonatal and paediatric patients mainly within clinical studies. The findings of these studies imply that lung monitoring using EIT provides valuable information on regional lung ventilation and aeration at the bedside that might improve the therapy and care of this fragile patient population. In view of this postulated future use of EIT in neonatology and paediatrics we have conducted an international survey to assess the perceived usefulness of several measures derived from EIT examinations. *Approach:* A questionnaire validating the clinical usefulness of 14 previously described EIT measures was designed and sent to 36 clinicians with previous experience with EIT in neonatal and paediatric patients. A numerical rating scale was used to assess the usefulness of each measure. *Main results:* 34 clinicians from 12 countries responded to the invitation and 32 filled in the questionnaire. The mean clinical and EIT experience ( $\pm$ SD) of the respondents was 19.4 $\pm$ 9.1 yr and 7.7 $\pm$ 5.8 yr, respectively. The top-rated measures were the global inhomogeneity index, silent spaces, change in end-expiratory lung impedance and ventrodorsal centre of ventilation. The bottom-rated were the regional respiratory time constant, tidal volume normalised to ml, respiratory rate and heart rate on the last rank. *Significance:* The survey revealed that EIT measures characterizing the ventilation and aeration distribution and the degree of their heterogeneity were deemed particularly useful. Respiratory rate, heart rate and overall tidal volume were considered less useful probably because these parameters are already routinely assessed by other conventional methods.

Keywords: EIT, electrical bioimpedance, respiratory system, lung, lung imaging

## 1. Introduction

Electrical impedance tomography (EIT) seems to be predestined for the clinical use in neonates, infants and children for the following reasons: EIT not only provides novel information on regional lung ventilation and aeration directly at the bedside that is not accessible to any other already established medical method but it also exhibits several features that make it particularly suitable for the use in patients at a young age. These features are mainly its non-invasiveness, radiation-free imaging and lack of need of patient's cooperation.

Nonetheless, the use of EIT in neonatal and paediatric patients is still limited to only a few centres in the world. EIT has been applied there in the delivery rooms, neonatal and paediatric intensive care units, operation theatres and pulmonary function labs with the aim of monitoring regional lung function and examining its disease- and treatment-related changes. EIT has been

used in neonatal and paediatric patients typically within clinical studies and not during routine care. The results of the published clinical studies show the potential of EIT chest monitoring in this patient population, as summarised in the recent consensus article (Frerichs *et al.*, 2017) and specifically highlighted in one of its online supplements dedicated to the EIT use in neonates, infants and children. The core patient population that might benefit from chest EIT are critically ill patients treated in neonatal (Chatziioannidis *et al.*, 2013; Miedema *et al.*, 2013; Rossi *et al.*, 2013; Hough *et al.*, 2014; van der Burg *et al.*, 2016) and paediatric intensive care units (Steinmann *et al.*, 2013; Wilsterman *et al.*, 2016; Lupton-Smith *et al.*, 2017; Dmytrowich *et al.*, 2018; Davies *et al.*, 2019). Children suffering from chronic lung diseases like cystic fibrosis or asthma (Lehmann *et al.*, 2016; Mueller *et al.*, 2018; Ngo *et al.*, 2018; Vogt *et al.*, 2018; Roethlisberger *et al.*, 2018), especially those in the pre-school age that cannot reliably perform conventional forced spirometric pulmonary function tests, represent another large patient group where EIT might become clinically relevant.

One of the reasons why the numbers of clinical studies in neonatology and paediatrics are relatively small compared with those in adult patients and why EIT is not routinely used in this setting yet was the unavailability of EIT devices approved for the use in neonates and infants by notified bodies. The introduction of devices certified for this specific patient population to the medical market and the availability of dedicated EIT electrode interfaces even for preterm neonates (Sophocleous *et al.*, 2018) can be expected to increase the clinical use. Another reason that still might be hampering the general use of EIT in the clinical setting is the lack of a consensus on the most adequate EIT measures needed for clinical decision-making. Although the recent consensus article on EIT (Frerichs *et al.*, 2017) provided a thorough summary and categorisation of EIT measures developed and used so far, it did not specify which parameters might be of particular relevance for the EIT use in neonatal and paediatric patients.

Therefore, we have conducted a survey among neonatologists and paediatricians active in EIT research with the aim of obtaining their personal assessment of multiple parameters that can be derived from EIT recordings.

## 2. Methods

### 2.1. Survey

Thirty six clinicians with variable degrees of experience with EIT from 13 countries were contacted by email. All of them have previously published articles on EIT use in neonates and children in peer-reviewed medical journals and/or participate in the currently running large clinical study on critically ill patients treated in neonatal and paediatric intensive care units (<http://cradlproject.org>). The authors of this article did not participate in the survey.


The contacted clinicians were asked to complete a questionnaire for evaluation of different EIT measures. The following 14 established EIT parameters were assessed: relative tidal volume ( $V_{T\text{rel}}$ ), absolute tidal volume normalised to ml ( $V_{T\text{abs}}$ ), ventrodorsal centre of ventilation ( $\text{CoV}_{\text{vd}}$ ) and ventilation ratio (v/d), right-to-left centre of ventilation ( $\text{CoV}_{\text{rl}}$ ) and ventilation ratio (r/l), change in end-expiratory lung impedance ( $\Delta\text{EELI}$ ), respiratory rate (RR), heart rate (HR), low tidal variation regions, i.e. so-called silent spaces (SS), regional respiratory system compliance ( $C_{\text{rs reg}}$ ), global inhomogeneity index (GI), coefficient of variation (CV) and regional respiratory time constants ( $\tau_{\text{reg}}$ ). The questionnaire included a brief description of each parameter, based on the definitions provided in the consensus article on chest EIT (Frerichs *et al.*, 2017), and its possible clinical applications. Clinicians were asked to evaluate the EIT parameters on a numerical rating scale ranging from 0 (not useful) to 7 (very useful).

The questionnaire additionally collected the information on the primary clinical specialty of the respondents (neonatology or paediatrics) and on the overall duration of both clinical and previous EIT experience in years. Subgroups of less ( $\leq 5$  years) and more experienced ( $> 5$  years) clinicians and EIT experts were then formed.

The participants in the survey were prompted to provide a short comment on each of the assessed EIT parameters in an allocated space below the ranking scale. At the end of the questionnaire, they were also asked whether there was any other EIT parameter they would wish to obtain from EIT examinations in neonates and children. Finally, an option for providing 'other' comments was given.

## 2.2. Data analysis

Data are shown as mean values  $\pm$  standard deviations. Histograms of the counts of each score's numerical values were generated and the percentages of positive and highly positive ratings of each EIT measure calculated. The ratings of all EIT parameters were compared with the repeated measures ANOVA with Bonferroni's test for multiple comparisons. The effects of the primary clinical specialty and the length of EIT experience on the ratings of each EIT measure were assessed by the unpaired t test. (The possible effect of the length of clinical experience was not statistically tested because of a very low number of less experienced colleagues participating in the survey, see below.) All statistical analyses were performed using GraphPad Prism 5.01 (GraphPad Software Inc., San Diego, CA, USA).

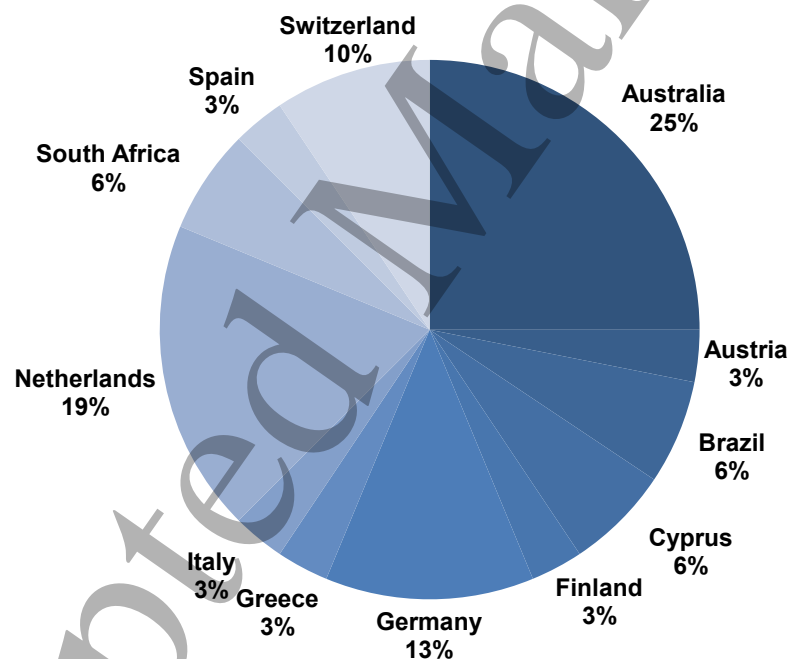

**Figure 1.** The countries of origin of 32 neonatologists and paediatricians who filled in the questionnaire on the clinical usefulness of EIT measures.


## 3. Results

### 3.1. Respondents

Thirty four clinicians responded to our initial invitation to fill in the questionnaires. Thirty two

colleagues with a mean clinical experience of  $19.4 \pm 9.1$  years and experience with EIT of  $7.7 \pm 5.8$  years from 12 countries on four continents completed the questionnaire (figure 1). Twenty respondents came from Europe, 8 from Australia, 2 from Africa and 2 from South America. Two colleagues declined the invitation because of illness and cessation of working as a medical doctor. Two of the contacted colleagues did not react in spite of repeated invitations.

Sixty six percent of the colleagues who filled in the questionnaires specified neonatology as their primary clinical specialty, the rest were paediatricians (figure 2, left). The majority of the clinicians participating in the survey (91 %) had an overall clinical experience of more than 5 years (figure 2, middle). EIT experience of more than 5 years was documented by 56 % of the survey participants (figure 2, right).

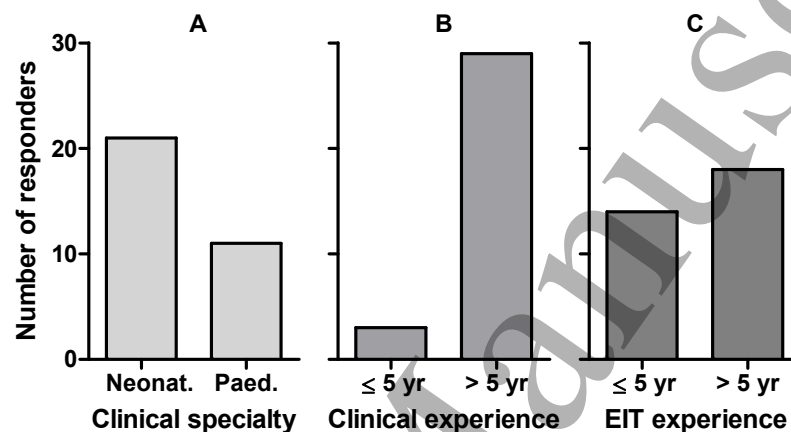

**Figure 2.** Numbers of respondents with the primary clinical specialty of neonatology (Neonat.) and paediatrics (Paed.) (A) with a short ( $\leq 5$  yr) and a long ( $> 5$  yr) clinical (B) and EIT experience (C).

The respondents had experience in the use of different EIT devices. Table 1 provides a list of EIT devices that were used by the survey participants in their previous EIT research. This list covers all devices applied in neonatal and paediatric EIT research in the past 20 years. 41 % of the respondents had experience with more than one device.

**Table 1.** List of EIT devices used by the survey participants

Manufacturer	EIT device	EIT interface	Number of EIT electrodes
CareFusion	Goe-MF II	individual electrodes	16
Dräger Medical GmbH	PulmoVista 500 <sup>†</sup>	electrode belt	16
Maltron International Ltd	Mark 1	individual electrodes	16
Maltron International Ltd	Mark 3.5 <sup>†</sup>	individual electrodes	8
SenTec AG *	BB <sup>2</sup> <sup>†</sup>	electrode belt	32
Timpel Medical	Enlight 1800 <sup>†</sup>	electrode stripes	32


<sup>†</sup> Commercially available as of March 2019, \* former Swisstom AG

### 3.2. Ratings

The ratings of the assessed EIT measures are summarised in figure 3. The scores differed significantly among the parameters ( $p = 0.0001$ ) with the highest ones given to GI ( $6.28 \pm 0.87$ ), SS ( $6.16 \pm 0.94$ ),  $\Delta EELI$  ( $6.16 \pm 0.87$ ) and  $CoV_{vd}$  ( $6.13 \pm 1.45$ ). The lowest rating was obtained for HR ( $4.19 \pm 2.26$ ,  $p < 0.01$  vs. all rated parameters except for  $V_{T\text{abs}}$ ). The other relatively low-rated EIT measures were RR ( $5.22 \pm 1.93$ ),  $V_{T\text{abs}}$  ( $5.34 \pm 1.81$ ) and  $\tau_{\text{reg}}$  ( $5.56 \pm 1.56$ ). There were no other statistically significant differences among the parameter ratings.

The clinical usefulness of all EIT parameters was consistently assessed as positive or even highly positive. The percentages of positive score values ranged from 69 % (HR) to 100 % ( $\Delta EELI$  and SS) (figure 3). The percentages of highly positive score values were in the range of 31 % (HR) to 88% (GI).

The primary clinical specialty exhibited no effect on the ratings of the majority of the assessed EIT parameters. Only  $\tau_{\text{reg}}$  was rated significantly more useful by the paediatricians ( $p = 0.0185$ ) and CV showed a trend for higher scores ( $p = 0.0655$ ) in the same group of specialists. The ratings were similarly not affected by the duration of EIT experience, except for v/d that was rated less positive by the clinicians with less years of EIT experience ( $p = 0.0283$ ) and  $V_{T\text{abs}}$  with a trend for worse ratings by the survey participants with more than 5 years of EIT experience ( $p = 0.0750$ ).

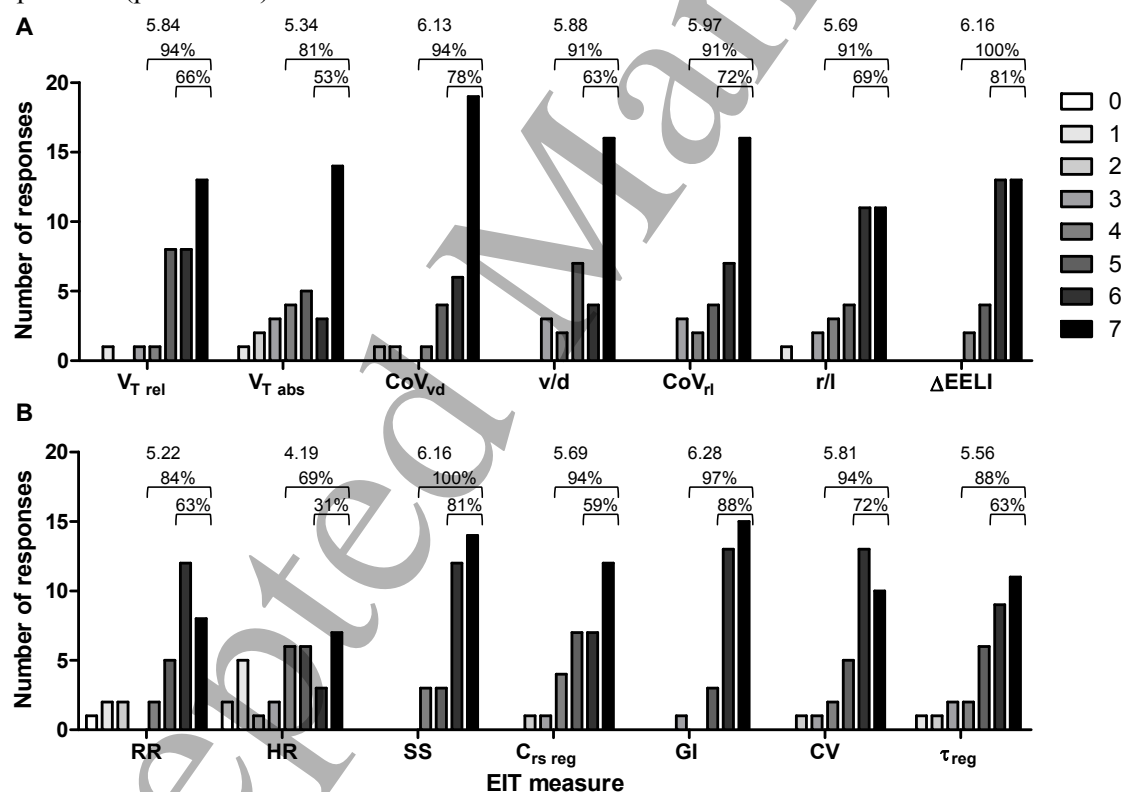

**Figure 3.** Assessment of clinical usefulness of EIT measures using a rating scale in the range of 0 (not useful, white columns) to 7 (very useful, black columns). The darker columns to the right of the ticks on the X axis imply clinical usefulness of the respective measure. The numbers above the columns in each panel

show the average values of the ratings (top row) and the percentages of positive (middle row) and highly positive responses (bottom row).  $V_{T\text{ rel}}$ , tidal volume determined by EIT in arbitrary units;  $V_{T\text{ abs}}$ , tidal volume determined by EIT normalised to global tidal volume in ml;  $\text{CoV}_{\text{vd}}$ , centre of ventilation in the ventrodorsal direction;  $v/d$ , ratio of ventrodorsal ventilation distribution;  $\text{CoV}_{\text{rl}}$ , centre of ventilation in the right-to-left direction;  $r/l$ , ratio of right-to-left ventilation distribution;  $\Delta\text{EELI}$ , change in end-expiratory lung impedance (A), RR, respiratory rate determined by EIT; HR, heart rate determined by EIT; SS, silent spaces, i.e. regions with low tidal impedance variation;  $C_{\text{rs reg}}$ , regional respiratory system compliance; GI, global inhomogeneity index; CV, coefficient of variation;  $\tau_{\text{reg}}$ , regional respiratory time constant (B).

### 3.3. Comments

Thirteen out of the 32 participants in the survey provided additional hand-written or typed comments on the assessed EIT measures, specified which additional EIT parameters they would prefer to obtain and/or wrote some 'other' comments.

Four to six comments were provided for each EIT parameter. The comments regarding  $V_{T\text{ rel}}$  and  $V_{T\text{ abs}}$  pointed out that conventional  $V_T$  measurements in neonates ventilated invasively through uncuffed tubes were unreliable and that EIT might become a reliable alternative. One respondent considered  $V_{T\text{ rel}}$  to be "one of the most important EIT measures" whereas another one wrote that " $V_{T\text{ rel}}$  change might be interesting as a trend but absolute reliable values are more interesting for real time titration". One participant in the survey wrote that " $V_{T\text{ abs}}$  would be most helpful for clinicians" but that he would be "extremely cautious to use it in clinical practice". He continued that "one had to be careful to 'calibrate' a chest slice as a representation of the whole lung" and that movement of a baby or changed body position would change the impedance characteristics and impact the value of  $V_{T\text{ abs}}$ . This is in harmony with the comment of another respondent who mentioned that " $V_{T\text{ abs}}$  would only be useful if there was sufficient data to confirm accuracy over a period of time (e.g. 4-6 hours)". One participant in the survey raised the point of "reliability of both  $V_{T\text{ rel}}$  and  $V_{T\text{ abs}}$  assessment obtained from EIT examinations with electrodes placed in only one plane". The clinical application of  $V_{T\text{ abs}}$  was considered "doubtful" by one respondent "because (transcutaneous) blood gases, like arterial partial pressure of  $\text{CO}_2$ , were more useful to monitor under/over-ventilation".

The measures of ventrodorsal ( $\text{CoV}_{\text{vd}}$ ,  $v/d$ ) and right-to-left ventilation distribution ( $\text{CoV}_{\text{rl}}$ ,  $r/l$ ) received overall positive comments.  $\text{CoV}_{\text{vd}}$  and  $v/d$  were considered valuable mainly for the assessment of atelectasis, recruitment and positive end-expiratory pressure (PEEP) titration. One respondent mentioned that these measures were "easily understood by clinicians" and that both CoV measures "were critical for commercial EIT as they provided a simple trend measure of regional  $V_T$  distribution and also allowed for identification of adverse events (e.g. endotracheal tube malposition and airleak)". The ability of  $\text{CoV}_{\text{rl}}$  and  $r/l$  to confirm suspected pneumothorax was described as "extremely useful" by another participant in the survey. One clinician added that more harmful chest X-ray was presently used to establish the pneumothorax diagnosis. Three colleagues commented that  $\text{CoV}_{\text{vd}}$  and  $v/d$  were similar parameters describing one phenomenon. Consequently, they recommended that only one of the two would be needed in the future. ( $\text{CoV}_{\text{rl}}$  and  $r/l$  received the same comment.) One participant in our survey wrote that "the larger the child, the more added value these parameters might deliver".

Four comments were given on  $\Delta\text{EELI}$ . Two respondents were convinced of the clinical relevance of this measure describing it as "very useful" and important for "guidance of recruitment manoeuvres at the bedside." However, one participant in the survey wrote that he could "not figure out how much this index would help during ventilation" and another wrote that "intrathoracic fluid volume changes might make it impossible to quantify end-expiratory lung volumes" using this measure.

1  
2  
3  
4  
5  
6  
7  
8  
9  
10  
11  
12  
13  
14  
15  
16  
17  
18  
19  
20  
21  
22  
23  
24  
25  
26  
27  
28  
29  
30  
31  
32  
33  
34  
35  
36  
37  
38  
39  
40  
41  
42  
43  
44  
45  
46  
47  
48  
49  
50  
51  
52  
53  
54  
55  
56  
57  
58  
59  
60

Regarding HR and RR it was stated that “ECG monitoring is standard” and that “other techniques like respiratory inductive plethysmography, electromyography and measurement of electrical chest trans-impedance through ECG electrodes are competitive to EIT”. Another participant in the survey pointed out that “whilst respiratory rate was measured accurately already in neonatal intensive care unit (NICU) using other systems, the ability to have less monitoring attached to a patient’s skin was desirable (and increased commercial viability of EIT systems)”. Two of the participants highlighted the importance of RR for the classification and therapy of apnoea, another one for the monitoring of spontaneously breathing or non-invasively ventilated patients. One colleague wrote that it would be desirable to have “reliable RR and HR data in the first couple of minutes after birth”.

SS received comments from five participants in our survey. One of them wrote that he expected SS to “guide clinicians during mechanical ventilation“ and show “if certain lung areas needed to be recruited or if a pneumothorax or chylothorax were present and if new air or fluid was accumulating”. Another noted that SS “might be useful since most other techniques were unable to detect these ‘silent zones’ of ventilation” although EIT might not cover the whole lung. One clinician considered SS to be useful “to confirm if the aims of used lung protective ventilation were fulfilled”. Although SS obtained only high scores, one comment implied that this parameter might need further refinement. This participant in the survey wrote that “from his experience, the clinicians understood this measure and found it very simple to visualize”. But he continued that “from his experimental and clinical observations he was not sure it was actually showing the clinically relevant information of whether a part of the lung was truly not engaging in much tidal change”. The comment further implied that the availability of specific neonatal/paediatric reconstruction models might improve the value of SS.

The comments on  $C_{rs\ reg}$  were interesting in the way that they all gave a positive feedback at the beginning but afterwards all of them mentioned potential limitations. For instance one respondent wrote that  $C_{rs\ reg}$  “would surely revolutionise the way ventilation was optimised” but he warned of “possible pitfalls related to real time information and measurement accuracy”. Another participant wrote that  $C_{rs\ reg}$  was an excellent measure for research but “less important for clinicians since most neonatologists did not even know what to do with this value”. The next respondent’s standpoint was that  $C_{rs\ reg}$  would “provide a true point of difference over other monitoring systems. It would be useful even if just relative (e.g. one region-of-interest vs. another)”. Finally one respondent correctly wrote that changes in  $C_{rs\ reg}$  reflected changes in regional  $V_T$  if the pressure swings were identical in the EIT plane.

The two global measures of ventilation heterogeneity (GI and CV) received very similar comments. One representative comment stated that “these single-number measures were what clinicians really wanted at this stage for simple understanding of temporal changes in ventilation distribution. Whether this was GI or CoV did not really matter.” One respondent commented that “these parameters were interesting but that there was some redundancy between them and the other rated measures”. He continued that “having too many parameters might ‘pollute’ the interpretation of findings”.

$\tau_{reg}$  was described as an “interesting measure for patient care and research” and even as “one of the most useful parameters”. However, one respondent wrote that although he personally valued this measure “in reality, it was an ‘expert’ user parameter as most clinicians did not understand simple time constants”. One participant in our survey wrote that  $\tau_{reg}$  might be a useful parameter for avoiding intrinsic PEEP during mechanical ventilation.

The prompt to provide a comment on the need of any additional EIT parameters was followed by three participants in the survey, all of them paediatricians. (None of them provided



any other comment on the other 14 EIT measures specified in the questionnaire.) They requested EIT measures allowing 1) the assessment of overdistension and collapse after each change in ventilator settings, 2) the early detection of pneumothorax and pulmonary interstitial emphysema and 3) the estimation of pulmonary perfusion and cardiac output.

The possibility to add 'other' comments at the end of the questionnaire was taken by four respondents. One paediatrician reported her positive experience with CoV and GI in ARDS patients and another one wrote down his wish to obtain EIT measures at various chest levels, i.e. also from apical and basal lung regions. One neonatologist commented that "she would have probably had more comments if her experience with EIT had been longer". Finally, one neonatologist pointed out that "clinicians would ideally need guidance on how to interpret all the EIT parameters". He gave an example of an alarm warning of a possible wrong location of the endotracheal tube.

#### 4. Discussion

Our survey focused exclusively on the clinical use of EIT in neonatal and paediatric patients. The first clinical chest EIT studies in neonates and children were initiated in the mid-1990s, much later than those in adults. These were initially mainly case reports documenting the feasibility and safety of EIT examinations in neonatal and paediatric patients with the first unique findings in specific diseases and clinical situations (Taktak *et al.*, 1996; Frerichs *et al.*, 1999; Frerichs *et al.*, 2001). The latest studies have already been performed on larger groups of patients with more sophisticated study protocols and EIT data analyses (van der Burg *et al.*, 2015; Hough *et al.*, 2016; Lehmann *et al.*, 2016; van der Burg *et al.*, 2016; Bhatia *et al.*, 2017; Mazzoni *et al.*, 2017; Lupton-Smith *et al.*, 2017; Vogt *et al.*, 2018). These research activities are driven by the presumed benefit of EIT monitoring in neonates and children with potentially high clinical relevance. This was also documented by several reviews describing the EIT use in these patients (Pillow *et al.*, 2006; Frerichs *et al.*, 2008; Riedel and Frerichs, 2010; Reiterer *et al.*, 2015; Smallwood and Walsh, 2017; Frerichs *et al.*, 2017).

Several commercial and non-commercial EIT devices have been used in neonates and children up to now. EIT examinations were performed at different scan rates and with dissimilar electrode interfaces (single electrodes and electrode belts). The acquired data was analysed using the built-in graphical user interfaces and/or dedicated, study-specific software. The analyses were often carried out offline. A variety of EIT measures was calculated from the EIT recordings. However, as the EIT technology advances and its use in a true clinical setting beyond the clinical studies becomes more probable, there is a need to define which EIT measures are relevant for clinical decision-making. The first step in defining these essential measures is the feedback from the current users of EIT regarding the clinical usefulness of the measures used so far. To obtain these perceptions on EIT measures and their use in neonates and children with ratings, was the primary motivation for our survey.

The users of EIT in neonatology and paediatrics are a relatively small group of researchers/clinicians but they are a very active and dedicated group within the EIT community. This was also confirmed by the high response rate to our invitation to participate in the survey of 94.4%. The analysis of the filled-in questionnaires revealed that the participants valued mainly those EIT measures that assessed the distribution of ventilation and aeration and its heterogeneity. Therefore, the top ranked measures were GI, SS,  $\Delta EELI$  and  $CoV_{vd}$ . The rankings of measures like HR, RR and  $V_{T\ abs}$  that are already routinely derived from conventional monitoring were the lowest. However, as pointed out by one of the participants, if these measures were reliably

1  
2  
3  
4  
5  
6  
7 provided by EIT then duplicate assessments by other techniques might become unnecessary. This  
8 might reduce the numbers of electrodes and other interfaces applied on the skin which is a  
9 relevant issue in preterm neonates and at birth.

10 The value and utility of the measures of regional respiratory system mechanics ( $C_{rs\ reg}$  and  
11  $\tau_{reg}$ ) was acknowledged but many respondents were concerned that these parameters might be too  
12 'complicated' for the future clinical users without a research background. However, such  
13 supposedly 'difficult' measures need not be explicitly shown and communicated to the clinical  
14 personnel in the later routine EIT monitoring. For example, one of the participants in the survey  
15 commented that the knowledge of regional overdistension and collapse would be valuable  
16 feedback information for assessing the effects of changed ventilator settings. To determine the  
17 degree of relative overdistension and collapse during a PEEP trial, changes in  $C_{rs\ reg}$  are typically  
18 calculated from the EIT recordings (Costa *et al.*, 2009) but the underlying calculations are not  
19 necessarily presented to the users.

20  
21 Interestingly,  $\tau_{reg}$  received higher ratings by the paediatricians. This might be related to the  
22 fact that older children are often treated for obstructive lung diseases where the utility of this  
23 measure of respiratory system mechanics might be more intuitive. In adult patients, the use of  $\tau_{reg}$   
24 for visualising the spatial heterogeneity of regional expiration in chronic obstructive lung disease  
25 was recently documented (Karagiannidis *et al.*, 2018). The capacity of  $\tau_{reg}$  to identify  
26 bronchospasm and the effects of bronchodilator therapy and suctioning was demonstrated in one  
27 paediatric patient (de la Oliva *et al.*, 2017). CV as an overall measure of ventilation heterogeneity  
28 also trended to get higher scores by paediatricians which probably reflected the higher usage of  
29 this parameter in studies on regional lung function outside the NICU setting. Another difference  
30 between the perceptions of neonatologists and paediatricians was that only the participants with  
31 the latter specialty expressed their concern about EIT findings not being representative for the  
32 whole lung. This seems not to be the case in neonatal patients (van der Burg *et al.*, 2014).  
33 Therefore, only the paediatricians articulated their need of obtaining EIT information at several  
34 chest levels.

35 We have no clear explanation why the respondents with a short EIT experience rated v/d  
36 worse than those with a long one.  $CoV_{vd}$ , a measure with similar information content, was not  
37 assessed differently by the experienced and less experienced users of EIT. The trend to give  
38 higher scores to  $V_{T\ abs}$  by the less experienced users seems to be easier to explain. These users  
39 prefer units (millilitre or litre) they are familiar with. They start to perceive the potential pitfalls  
40 of calibrating the EIT signal with the global volumetric measurements only with increasing  
41 expertise in EIT. This assessment is supported by some of the rather cautious comments  
42 regarding this measure requesting that the long-term reliability of  $V_{T\ abs}$  needed to be proved.

43 Several comments of the respondents emphasised that the usefulness and clinical acceptance  
44 of EIT measures would decisively depend on whether the users of EIT technology were provided  
45 with adequate guidance how to interpret the findings. Thus, the essential issue is the generation of  
46 interpretation schemes for different EIT parameters and their combinations. Ideally, they should  
47 also take other clinical findings into account. As stressed by one of the participants in the survey,  
48 the clinicians need to be alarmed and trained that a highly asymmetric ventilation distribution  
49 might be indicative of a pneumothorax. Whether the alarm is triggered by  $CoV_{rl}$  or  $r/l$  or by a  
50 combination with other measures (e.g.  $\Delta EELI$ ) will depend on the underlying interpretation  
51 algorithm.

52  
53 The comments of the participants in the survey also draw attention to the potential utility of  
54 some measures derived from the EIT recordings in detecting adverse events like endotracheal  
55 tube malposition and pneumothorax or other pathological findings with high ventilation  
56  
57  
58  
59  
60

1  
2  
3  
4  
5  
6  
7  
8  
9  
10  
11  
12  
13  
14  
15  
16  
17  
18  
19  
20  
21  
22  
23  
24  
25  
26  
27  
28  
29  
30  
31  
32  
33  
34  
35  
36  
37  
38  
39  
40  
41  
42  
43  
44  
45  
46  
47  
48  
49  
50  
51  
52  
53  
54  
55  
56  
57  
58  
59  
60

asymmetry. Several published case reports seem to confirm the validity of this assessment (Frerichs *et al.*, 2005; Miedema *et al.*, 2011; Tingay *et al.*, 2019). The use of EIT in these situations might lead to faster finding of the proper diagnosis, treatment initiation and/or adjustment.

Our study has a few limitations that we wish to acknowledge. Although we successfully reached almost all clinical users of EIT in neonatology and paediatrics, the overall number of participants was still rather low. Our survey focused merely on the EIT measures and no assessment of possible technical limitations of EIT use in neonatal and paediatric patients was surveyed. The individual ratings of the measures and the comments may have been influenced by the personal experience of the users with specific EIT technology. The users normally use the EIT devices of just one brand and are not familiar with other devices. This was the case of 59 % of respondents in the survey. This may have induced a possible bias toward measures offered by the device available in their department or by the local research software used for offline data analysis. Finally, to establish a set of clinically relevant EIT measures undeniably not only subjective clinicians' assessments but also evidence-based data will be needed in the future.

## 5. Conclusion

A group of neonatologists and paediatricians with various degrees of experience with EIT have rated the clinical usefulness of a selected set of 14 EIT measures using a questionnaire. The ratings revealed that the parameters characterising ventilation and aeration distribution (GI, SS,  $\Delta EELI$  and  $CoV_{v,d}$ ) were considered particularly useful for neonatal and paediatric patients. The free-text comments written by the participants in this survey provided additional valuable information on the advantages, drawbacks and possible future directions in neonatal and paediatric EIT research, clinical implementation and use.

## Acknowledgments

We would like to express our gratitude to all the colleagues who participated in the survey. We also acknowledge the funding by the European Union's Horizon 2020 Research and Innovation Programme (Project CRADL, grant agreement number 668259).

## References

- Bhatia R, Davis P G and Tingay D G 2017 Regional volume characteristics of the preterm infant receiving first intention continuous positive airway pressure *J Pediatr* **187** 80-8 e2
- Chatziioannidis I, Samaras T, Mitsiakos G, Karagianni P and Nikolaidis N 2013 Assessment of lung ventilation in infants with respiratory distress syndrome using electrical impedance tomography *Hippokratia* **17** 115-9
- Costa E L, Borges J B, Melo A, Suarez-Sipmann F, Toufen C, Jr., Bohm S H and Amato M B 2009 Bedside estimation of recruitable alveolar collapse and hyperdistension by electrical impedance tomography *Intensive Care Med.* **35** 1132-7
- Davies P, Yasin S, Gates S, Bird D and Silvestre C 2019 Clinical scenarios of the application of electrical impedance tomography in paediatric intensive care *Sci Rep* **9** 5362
- de la Oliva P, Waldmann A D, Bohm S H, Verdu-Sanchez C, Perez-Ferrer A and Alvarez-Rojas E 2017 Bedside breath-wise visualization of bronchospasm by electrical impedance tomography could improve perioperative patient safety: A Case Report *A A Case Rep* **8** 316-9

- 1  
2  
3  
4  
5  
6  
7 Dmytrowich J, Holt T, Schmid K and Hansen G 2018 Mechanical ventilation guided by electrical  
8 impedance tomography in pediatric acute respiratory distress syndrome *J Clin Monit*  
9 *Comput* **32** 503-7
- 10 Frerichs I, Amato M B, van Kaam A H, Tingay D G, Zhao Z, Grychtol B, Bodenstern M, Gagnon  
11 H, Bohm S H, Teschner E, Stenqvist O, Mauri T, Torsani V, Camporota L, Schibler A,  
12 Wolf G K, Gommers D, Leonhardt S and Adler A 2017 Chest electrical impedance  
13 tomography examination, data analysis, terminology, clinical use and recommendations:  
14 consensus statement of the TRanslational EIT developmeNt stuDY group *Thorax* **72** 83-  
15 93
- 16 Frerichs I, Hahn G, Schiffmann H, Berger C and Hellige G 1999 Monitoring regional lung  
17 ventilation by functional electrical impedance tomography during assisted ventilation  
18 *Ann N Y Acad Sci* **873** 493-505
- 19 Frerichs I, Schiffmann H, Hahn G, Dudykevych T, Just A and Hellige G 2005 Funktionelle  
20 elektrische Impedanztomographie - eine Methode zur bettsetigen Ueberwachung der  
21 regionalen Lungenfunktion. *Intensivmed* **42** 66-73
- 22 Frerichs I, Schiffmann H, Hahn G and Hellige G 2001 Non-invasive radiation-free monitoring of  
23 regional lung ventilation in critically ill infants *Intensive Care Med.* **27** 1385-94
- 24 Frerichs I, Scholz J and Weiler N 2008 Regional lung function in critically ill neonates: A new  
25 perspective for electrical impedance tomography. In *Yearbook of Intensive Care and*  
26 *Emergency Medicine* Vincent J-L ed Springer-Verlag: Berlin pp 224-34
- 27 Hough J, Trojman A and Schibler A 2016 Effect of time and body position on ventilation in  
28 premature infants *Pediatr Res* **80** 499-504
- 29 Hough J L, Shearman A D, Liley H, Grant C A and Schibler A 2014 Lung recruitment and  
30 endotracheal suction in ventilated preterm infants measured with electrical impedance  
31 tomography *J Paediatr Child Health* **50** 884-9
- 32 Karagiannidis C, Waldmann A D, Roka P L, Schreiber T, Strassmann S, Windisch W and Bohm  
33 S H 2018 Regional expiratory time constants in severe respiratory failure estimated by  
34 electrical impedance tomography: a feasibility study *Crit. Care* **22** 221
- 35 Lehmann S, Leonhardt S, Ngo C, Bergmann L, Ayed I, Schrading S and Tenbrock K 2016 Global  
36 and regional lung function in cystic fibrosis measured by electrical impedance  
37 tomography *Pediatr Pulmonol* **51** 1191-9
- 38 Lupton-Smith A, Argent A, Rimensberger P, Frerichs I and Morrow B 2017 Prone positioning  
39 improves ventilation homogeneity in children with acute respiratory distress syndrome  
40 *Pediatr Crit Care Med* **18** e229-e34
- 41 Mazzoni M B, Perri A, Plebani A M, Ferrari S, Amelio G, Rocchi A, Consonni D, Milani G P  
42 and Fossali E F 2017 Electrical impedance tomography in children with community  
43 acquired pneumonia: preliminary data *Respir Med* **130** 9-12
- 44 Miedema M, Frerichs I, de Jongh F H, van Veenendaal M B and van Kaam A H 2011  
45 Pneumothorax in a preterm infant monitored by electrical impedance tomography: a case  
46 report *Neonatology* **99** 10-3
- 47 Miedema M, van der Burg P S, Beuger S, de Jongh F H, Frerichs I and van Kaam A H 2013  
48 Effect of nasal continuous and biphasic positive airway pressure on lung volume in  
49 preterm infants *J Pediatr* **162** 691-7
- 50 Mueller J L, Muller P, Mellenthin M, Murthy R, Capps M, Alsaker M, Deterding R, Sagel S D  
51 and DeBoer E 2018 Estimating regions of air trapping from electrical impedance  
52 tomography data *Physiol. Meas.* **39** 05NT1
- 53  
54  
55  
56  
57  
58  
59  
60

- 1  
2  
3  
4  
5  
6  
7 Ngo C, Dippel F, Tenbrock K, Leonhardt S and Lehmann S 2018 Flow-volume loops measured  
8 with electrical impedance tomography in pediatric patients with asthma *Pediatr Pulmonol*  
9 **53** 636-44
- 10 Pillow J J, Frerichs I and Stocks J 2006 Lung function tests in neonates and infants with chronic  
11 lung disease: global and regional ventilation inhomogeneity *Pediatr Pulmonol* **41** 105-21
- 12 Reiterer F, Sivieri E and Abbasi S 2015 Evaluation of bedside pulmonary function in the neonate:  
13 From the past to the future *Pediatr Pulmonol* **50** 1039-50
- 14 Riedel T and Frerichs I 2010 Electrical impedance tomography *Eur Respir Mon* **47** 195-205
- 15 Roethlisberger K, Nyilas S, Riedel T, Wolfensberger J, Singer F and Latzin P 2018 Short-term  
16 effects of elastic chest wall restriction on pulmonary function in children with cystic  
17 fibrosis *Respiration* **96** 535-42
- 18 Rossi F d S, Yagui A C, Haddad L B, Deutsch A D and Rebello C M 2013 Electrical impedance  
19 tomography to evaluate air distribution prior to extubation in very-low-birth-weight  
20 infants: a feasibility study *Clinics (Sao Paulo)* **68** 345-50
- 21 Smallwood C D and Walsh B K 2017 Noninvasive monitoring of oxygen and ventilation *Respir*  
22 *Care* **62** 751-64
- 23 Sophocleous L, Frerichs I, Miedema M, Kallio M, Papadouri T, Karaoli C, Becher T, Tingay D  
24 G, van Kaam A H, Bayford R and Waldmann A D 2018 Clinical performance of a novel  
25 textile interface for neonatal chest electrical impedance tomography *Physiol. Meas.* **39**  
26 044004
- 27 Steinmann D, Engehausen M, Stiller B and Guttmann J 2013 Electrical impedance tomography  
28 for verification of correct endotracheal tube placement in paediatric patients: a feasibility  
29 study *Acta Anaesthesiol. Scand.* **57** 881-7
- 30 Taktak A, Spencer A, Record P, Gadd R and Rolfe P 1996 Feasibility of neonatal lung imaging  
31 using electrical impedance tomography *Early Hum Dev* **44** 131-8
- 32 Tingay D G, Waldmann A D, Frerichs I, Ranganathan S and Adler A 2019 Electrical impedance  
33 tomography can identify ventilation and perfusion defects: a neonatal case *Am. J. Respir.*  
34 *Crit. Care Med.* **199** 384-6
- 35 van der Burg P S, de Jongh F H, Miedema M, Frerichs I and van Kaam A H 2016 Effect of  
36 minimally invasive surfactant therapy on lung volume and ventilation in preterm infants *J*  
37 *Pediatr* **170** 67-72
- 38 van der Burg P S, Miedema M, de Jongh F H, Frerichs I and van Kaam A H 2014 Cross-sectional  
39 changes in lung volume measured by electrical impedance tomography are representative  
40 for the whole lung in ventilated preterm infants *Crit. Care Med.* **42** 1524-30
- 41 van der Burg P S, Miedema M, de Jongh F H, Frerichs I and van Kaam A H 2015 Changes in  
42 lung volume and ventilation following transition from invasive to noninvasive respiratory  
43 support and prone positioning in preterm infants *Pediatr Res* **77** 484-8
- 44 Vogt B, Lohr S, Zhao Z, Falkenberg C, Ankermann T, Weiler N and Frerichs I 2018 Regional  
45 lung function testing in children using electrical impedance tomography *Pediatr*  
46 *Pulmonol* **53** 293-301
- 47 Wilsterman M E, de Jager P, Blokpoel R, Frerichs I, Dijkstra S K, Albers M J, Burgerhof J G,  
48 Markhorst D G and Kneyber M C 2016 Short-term effects of neuromuscular blockade on  
49 global and regional lung mechanics, oxygenation and ventilation in pediatric acute  
50 hypoxemic respiratory failure *Ann Intensive Care* **6** 103
- 51  
52  
53  
54  
55  
56  
57  
58  
59  
60